



PRODUCT BROCHURE · HAIER BIOMEDICAL

Haier DW-25L262 -25°C Biomedical Lab Freezer , 115V

SKU DW-25L262 · Cold Storage · Low Temp Freezers



QUICK VIEW

The Haier DW-25L262 is an upright -25°C biomedical laboratory freezer with 262 L (9.3 cu ft) of storage, designed for the everyday -20°C-class needs

SOURCING

**Authorized U.S.
distributor**

Distributed by Legacy Lab Supplies

legacylabsupplies.com · (855) 746-7557 / 855-SHOP-LLS · quotes in 24 hours






Legacy Lab Supplies

Haier Biomedical
Intelligent Protection of Life Science

-25°C Biomedical Freezer

— DW-25L262

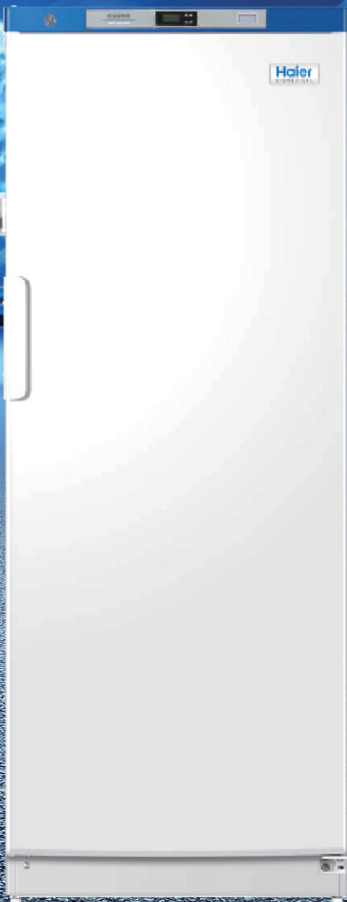

Hydrocarbon (HC) Energy-saving

DW-25L262 adopts an efficient hydrocarbon (HC) compressor and optimized refrigeration system, with power consumption as low as 0.73KW·h/24h.



Environmentally Friendly

The HC system is a completely green refrigeration solution. With LBA foamed polyurethane insulation in line with European ROHS directive for environmentally friendly materials.



Low Noise

Optimized hydrocarbon system and structural noise reduction design, DW-25L262 noise 34 dB, super quiet and low noise.



Security Lock

Mechanical lock plus lock catch, the mechanical door lock can be opened with one key, and the door lock can be padlocked. Secure and safe, suitable for multiple users.

Qingdao Haier Biomedical Co.,Ltd.

No.280 Feng Yuan Road, High-tech Zone,
Qingdao, 266109, P.R. China
Tel: +86-0532-88935593
Website: www.haiermedical.com



Haier Biomedical
International



Haier Biomedical
International



@haiermedicalint



@haierbiomedicalint



Haier Biomedical
International

-25°C Biomedical Freezer DW-25L262

Legacy Lab Supplies




Product Features



Microprocessor

- Microprocessor electronic temperature controller, LCD digital temperature display, display precision is 0.1° C, interior temperature set point -10~-25°C.
- Multiple alarms including high and low temperature, sensor error, power failure, door ajar, high ambient and remote alarm.
- Alarm mode: sound and light alarm.
- It can be connected to remote alarm.

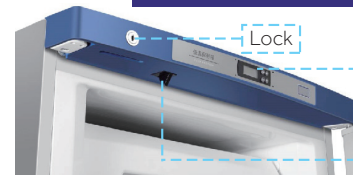

Lock

The new thermostat displays precision is 0.1

Door switch



Hydrocarbon System

Increasing cooling rate, reducing energy consumption and lowering noise



Porthole

There is a standard access port (25mm inside diameter) on the left side of the unit, which is convenient for users to test independently.



Multiple Interface Options

- Standard configuration: Remote alarm interface.
- Optional configuration: RS485.



Lock catch

Equipped with remote alarm terminal



Key W0=TOTAL WIDTH W1=FREEZER WIDTH W2=INNER CABINET WIDTH W3=INNER CABINET SHELF WIDTH D0=TOTAL DEPTH D1=INNER CABINET DEPTH D2=INNER CABINET SHELF DEPTH H0=TOTAL HEIGHT H1=INNER CABINET HEIGHT									
CODE MODEL	W0	W1	W2	W3	D0	D1	D2	H0	H1
DW-25L262	700mm	657mm	480mm	450mm	705mm	465mm	313mm	1730mm	1430mm

Specifications

Defrost Mode	Refrigerant	Sound Level (dB(A))	Cooling Performance (°C)	Temperature Range (°C)	Power Supply (V/Hz)	Power (W)
Manual	HC	34	-25	-10~-25	115/60	135
Electrical Current (A)	Power Consumption (kW·h/24h)	Capacity (L/Cu.Ft)	Net/Gross Weight (approx)(kg)	Internal Dimension (W*D*H)(mm)		
2.6	0.73	262/9.3	88/93	480*465*1430		
Exterior Dimension (W*D*H)(mm)	Container Load (20'/40'/40'H)	Alarms			Certification	
700*705*1730	21/45/45	High/Low Temperature, Remote Alarm, Power Failure, Sensor Error, High Ambient Temperature, Door Ajar			UL	

Haier Biomedical